

Secondary Set, Drip Chamber | REF C60




Secondary Set, Drip Chamber

This non-vented Secondary Set is used for leak proof and air tight transfer of hazardous drugs. A flow stop facilitates priming of the Secondary Set from the bag prior to adding the drug. The universal spike, equipped with the PhaSeal Connector, enables safe transfer of drug to the IV bag.

BILL OF MATERIALS

Pos.	Components	Material	Characteristics	Qty
1	Spike	ABS	White	1
2	Connector	ABS	White	1
3	Membrane	Elastomer	White	1
4	Drip chamber	DEHP free PVC	Transparent	1
5	Tubing	DEHP free PVC	Transparent	1
-	Protective spike cap	Polyethylene	Semi-transparent	1
-	Roller clamp	ABS	Blue	1
-	Male luer-lock	ABS	Transparent	1
-	Protective cap with Flow-stop	Polyethylene	Semi-transparent	1

ASSEMBLY DRAWING

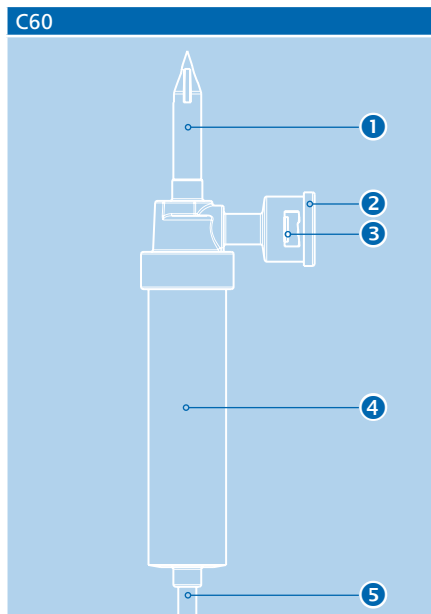

PHYSICAL CHARACTERISTICS

REF / Art. no.	Dimension length	Priming volume	Connector internal volume	Drip chamber, drops per volume	Packaging: User box	Packaging: Transport box
C60	102 cm (40 in.)	12.5 ml	0.15 ml	20 drops/ml	20 pcs	80 pcs